


## Product Information


# LabMaster – module

Video recording

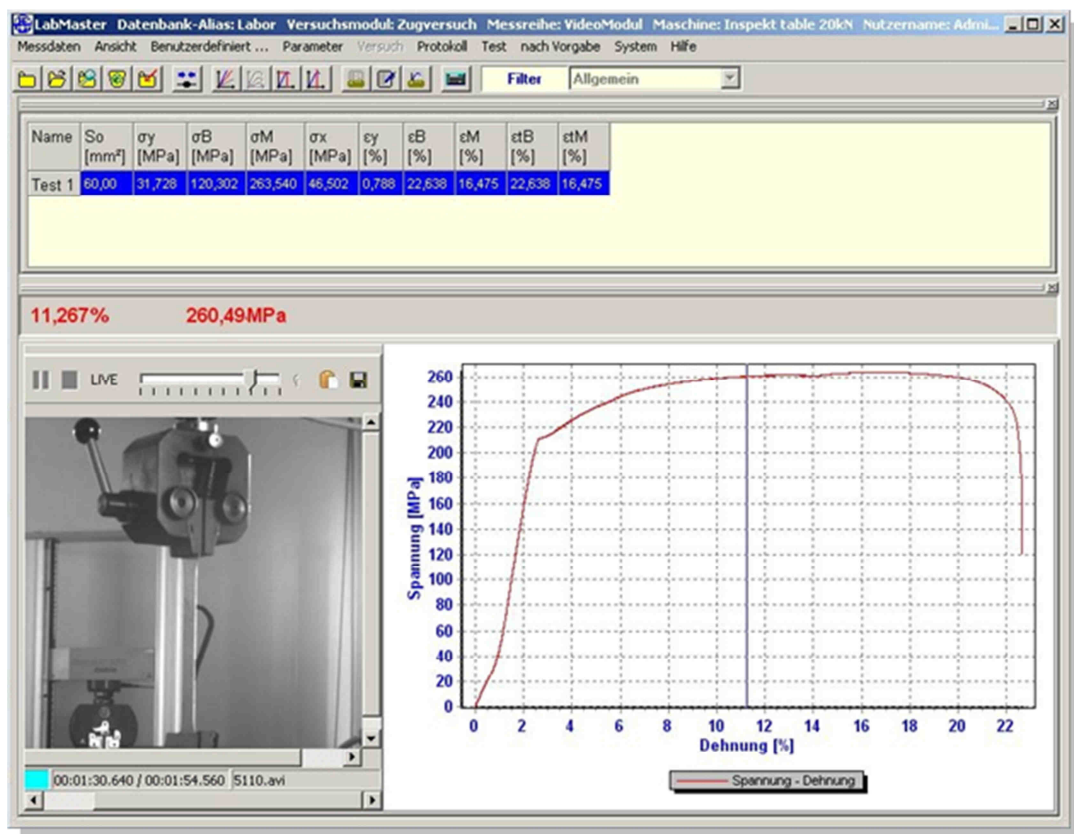

The module video recording allows the recording, storing and play back of video sequences of static or quasi static tests.

Loads and positions measured by the testing machine can be assigned to the video sequences time synchronously. Then, video stills can be displayed by moving the mouse cursor over the diagram. The image can be copied to the clipboard or saved to the hard disk.

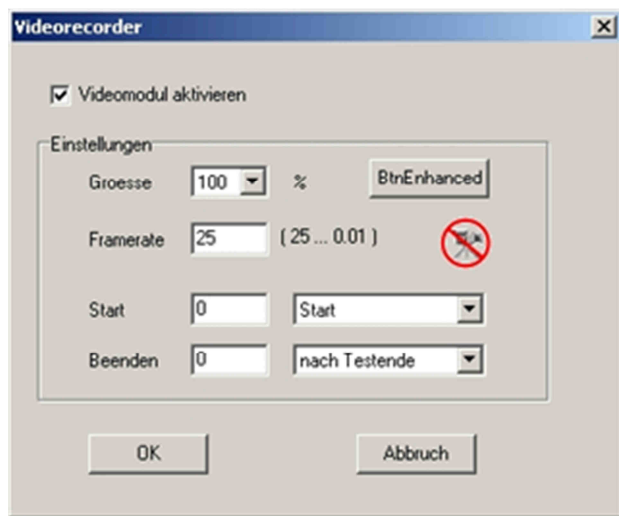

The simultaneous operation of the module video recording and of the testing software can only be done by a high-capacity computer provided by Hegewald & Peschke MPT GmbH. There are two options available for this module: either without a camera, which is interesting if the customer has a camera already. However, it must be ensured beforehand that the camera is suitable by providing it to Hegewald & Peschke MPT GmbH. The second option includes the digital camera and accessories (lenses, tripod).

- 18-014-012  
Material testing software LabMaster - module video recording without accessories (camera, lenses and tripod are provided by the customer)
- 18-014-013  
Material testing software LabMaster - module video recording incl. accessories (camera, lenses and tripod)